

January 2024


TECHNICAL PRODUCT BULLETIN TB-2024-01-01 Rev A

Subject: Shipper Core Changes


During the year of 2024, IC Biomedical will be making a design change to the shipper cores in all 10DX and Dual units. The previous design uses an aluminum mesh screen. The new design uses a rolled aluminum sheet with larger holes to allow for LN2 absorption. Functionally, there is no difference in absorption rate, NER, or hold time between the two designs.

Shipper Core Designs:

Previous Design

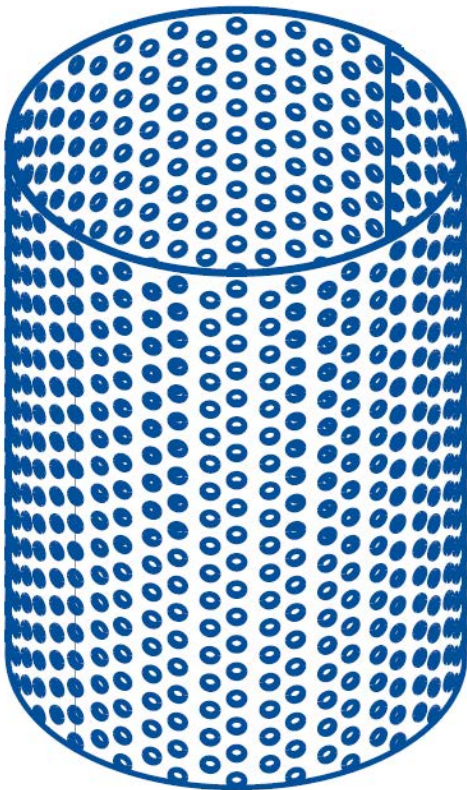

New Design

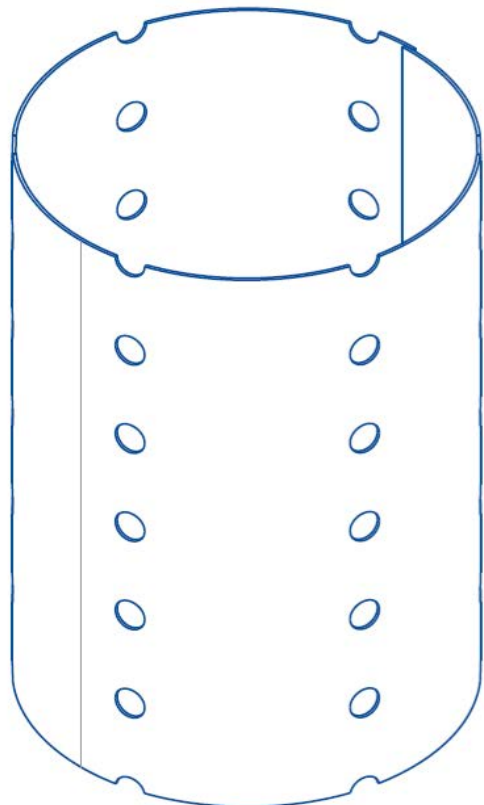

 1-855-750-8191

 www.icbiomedical.com

 1031 Cass White Road, NW
Cartersville, Georgia 30121